

# AUTOMATIC CAPSULE SEPARATING



## DESCRIPTION

This is the world first machine designed to separate capsules and re-use powder, granules, pellets, or micro-tablets which is built to replace the traditional way of milling, sieving and retrieving powder.


Our machine is first of its kind to separate and retrieve round pellets and micro-tablets. It is ergonomically designed; suitable for all capsule sizes with applicable change parts; easy to clean and maintain; user-friendly and meet the cGMP Standard in the industry.

It is excellent machine to own by any capsule manufacturing plants or R&D laboratories, where powders are extremely expensive and cGMP is strictly required.

## Principle

The capsules enter through the guide from the sorting system are pushed into separating device by the horizontal fork. The blades on the device hold the cap and body of the moving capsules.

The capsules moving according to the opening rail formed by the blades; then, it will separate the body and cap of capsule completely.

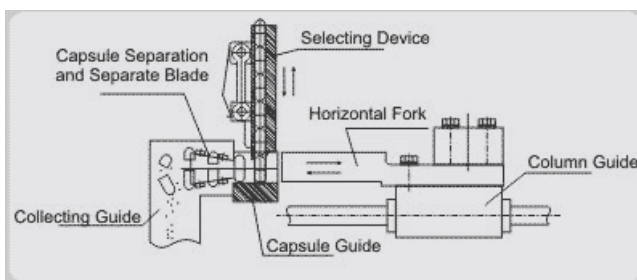

## DATA TECHNICAL

|                 |   |
|-----------------|---|
| Capsules Sizes: | 0# 1# 2# 3# 4#                            |
| Output          | 18,000 pcs/hour                           |
| Power           | ~220V 50Hz/60Hz 1P<br>~ 220V 50Hz/60Hz 3P |
| Dimension       | 60X600X1100mm                             |
| Weight          | 40 Kg                                     |

Your Local Distributor :

**JOTAM**

Telp: 021 4513089

Fax : 021 451 5743

Email: jotam@cbn.net.id